


VIPER® Cortical Fix Fenestrated Screw System

Surgical Technique
for OPEN & MIS Procedures

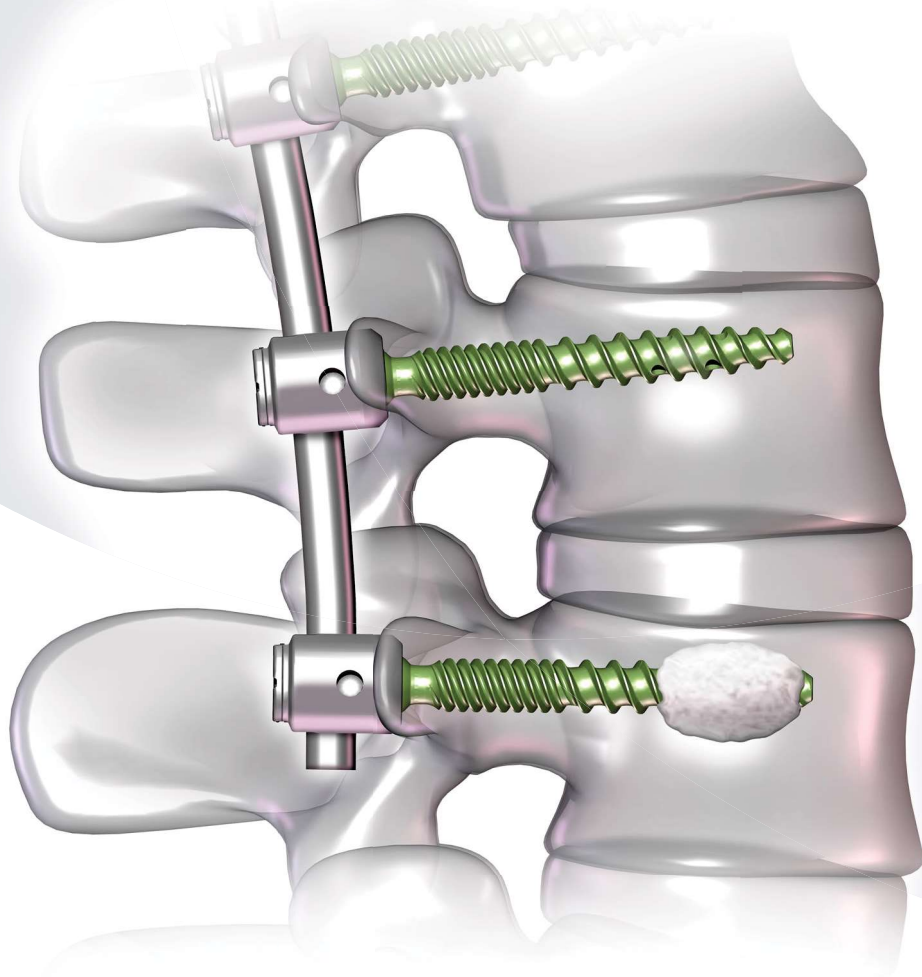

Features and Benefits

TOP NOTCH™ Feature

The TOP NOTCH Feature enables instruments from the EXPEDIUM System and the VIPER System to connect directly to the implants and simplifies intraoperative maneuvers.

Constant Thread Lead

The innovative thread form allows for a smooth transition between the dual and quad lead thread during insertion.

Cortical Thread Form

Cortical Thread Form provides enhanced pedicle fixation.

Fenestrations

1.6 – 1.75 mm fenestration diameters allow for optimal cement delivery and screw fixation.

Optimized Tip

Self-tapping and self-centering tip allows for insertion with or without tapping.

Cannulation



1.75 mm diameter allows for injection of cement through the screw.

TOP NOTCH Feature

Cortical Thread Form

Constant Thread Lead



Optimized Tip

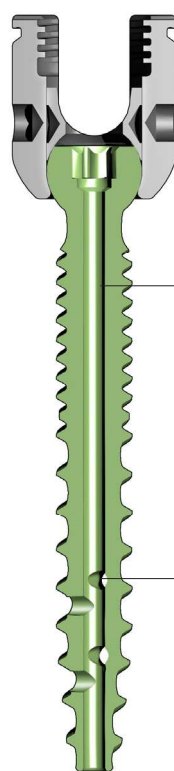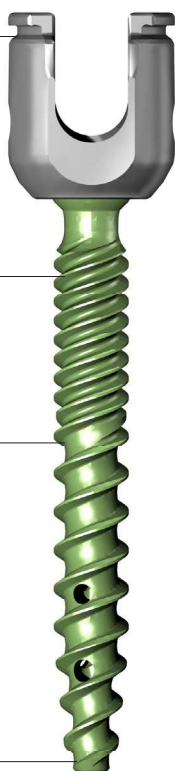



Cannulation

Fenestrations



Measurement Guide

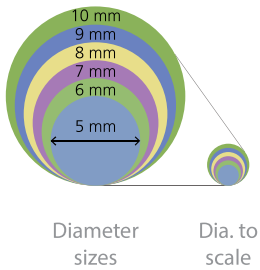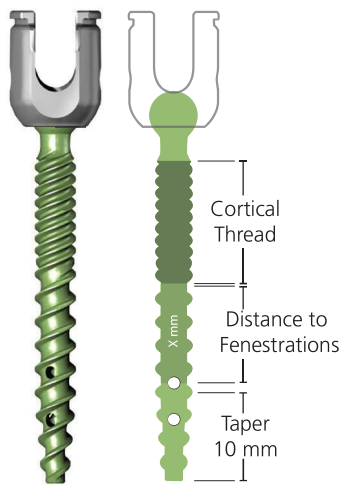

3 Fenestrations

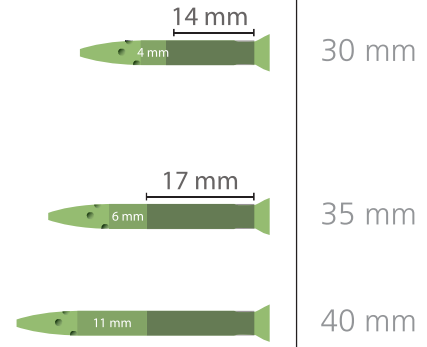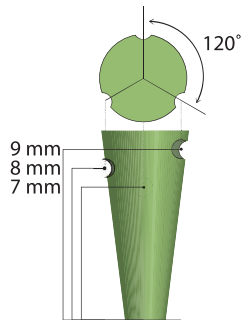

6 Fenestrations

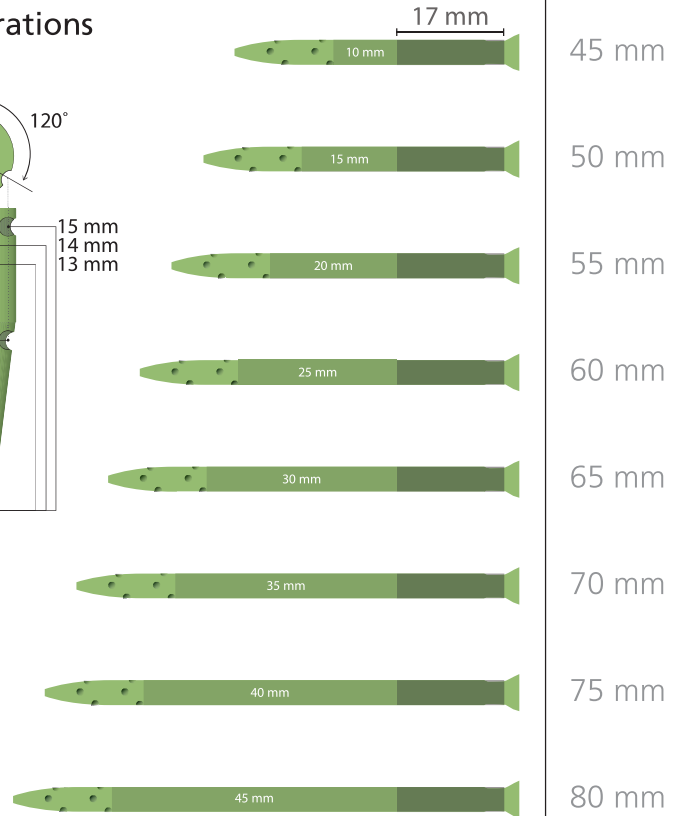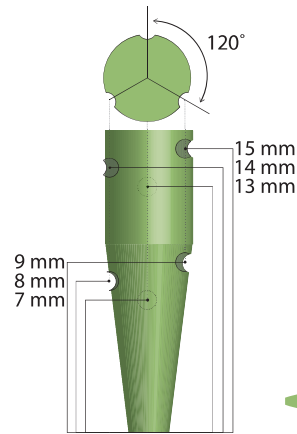

All values are rounded to the nearest whole number.